

OLYMPUS®



EVIS EXERA COLONOVideoscope

OLYMPUS CF TYPE 2T160L/I

POWERFUL DUAL-CHANNEL PERFORMANCE

Featuring two independent channels with diameters of 3.7 mm and 3.2 mm, this scope offers the most powerful treatment and suction capability for colonoscopic procedures




Main Features for Powerful Performance

- Dual channel design, with a 3.7 mm diameter channel for suction, and a 3.2 mm diameter channel for insertion of an Endo-Therapy accessory, ensures powerful suction capability and allows both channels to be used for Endo-Therapy accessories.
- 3 mm minimum depth of field for close-up observation.
- Extra-wide 140° field of view allows easier observation of a wider area of the site.
- 4-way angulation (180° up/down and 160° right/left) enables complete and comprehensive examination of the colon.
- 13.7 mm outer diameter insertion tube.
- Compatible with laser probes for expanded treatment versatility.
- Fully compatible with the CV-160/140/100.
- Scope ID function stores individual scope information and displays it on the monitor to facilitate endoscopy suite management.



Specifications

Optical System	Field of view	140°
	Direction of view	0° Forward viewing
	Depth of field	3 to 100 mm
Distal End	Outer diameter	13.8 mm
	Insertion Tube	Outer diameter
Bending Section	Angulation range	Up 180°, Down 180°, Right 160°, Left 160°
Working Length		L:1680 mm, l:1330 mm
Total Length		L:2005 mm, l:1655 mm
Instrument Channel	Inner diameter	A:3.7 mm, B:3.2 mm
	Minimum visible distance	5 mm from the distal end
	Endo-Therapy accessory entrance/exit position in field of view	

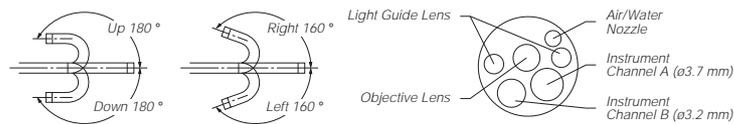

Olympus business areas

	Medical and health-care area
	Imaging and information area
	Industrial applications area

OLYMPUS®

FOCUS ON LIFE

OLYMPUS OPTICAL CO., LTD.
 2951 Ishikawa-cho, Hachioji-shi, Tokyo 192-8507, Japan
OLYMPUS OPTICAL CO. (EUROPA) GMBH.
 Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany
OLYMPUS AMERICA INC.
 2 Corporate Center Drive, Melville, N.Y. 11747-3157, U.S.A.
OLYMPUS KEYMED
 Keymed House, Stovok Road, Southend-on-Sea, Essex SS2 5QH, England
OLYMPUS SINGAPORE PTE LTD.
 491B River Valley Road #12-01/04 Valley Point Office Tower, Singapore 248373
OLYMPUS HONG KONG AND CHINA LIMITED.
 Room 1520-1527, Ocean Centre, 5 Canton Road, Tsimshatsui, Kowloon, Hong Kong
OLYMPUS BEIJING REPRESENTATIVE OFFICE
 Rm.818, South Tower, Beijing Kerry Centre, No.1 Guanghua Road, Chaoyang District, Beijing, 100020, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
 117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia
OLYMPUS AUSTRALIA PTY. LTD.
 1/104 Ferntree Gully Road, Oakleigh, VIC 3166, Australia

<http://www.olympus.com>